
REUDYN

USER MANUAL

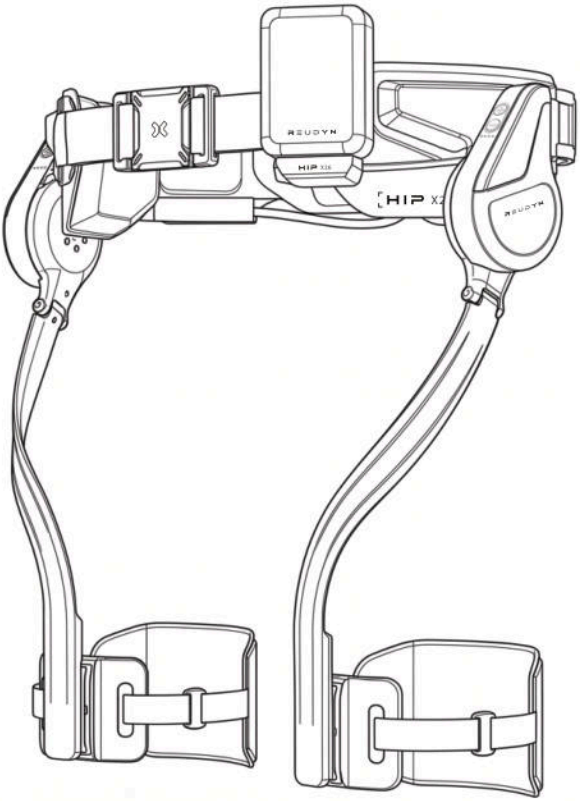CONTENTS

What's in the box---- 07
Product Overview---- 08
Diagram---- 08
Specifications---- 09
Wearing and Removal---- 09
Wearable Device---- 09
Removing the Device---- 14
Stowing the Device---- 15
First Use---- 17
First-Time Charging---- 17
Power On/Off---- 18
Connecting the Device---- 18
Smart Mode---- 19
Standard Mode---- 19
Custom Mode---- 19
Operating the Device---- 19
Control Buttons---- 19
Power On/Off---- 19
Switching Assistance Modes---- 19
Adjusting Assistance level---- 19
Rerunning wear detection---- 20
States Indicators---- 20
Battery---- 20
Removing the battery---- 21
Low Battery Safety System---- 22
Reudyn APP---- 22
Device Status---- 22
Firmware Updates---- 22
Maintenance and Storage---- 23
Cleaning the Main Unit---- 23
Battery Maintenance---- 24
Storage and Transportation---- 24
Storage---- 24
Transportation---- 24

⚠WARNING!

Indicates a hazardous situation that, if not avoided, could result in serious personal injury.

Do not use this product if you are pregnant; have skin lesions on areas that contact the product; have an infectious disease; are recovering from surgery; have a fracture; have undergone coronary artery bypass surgery; or have internal metal implants.

Before inserting the charging cable, ensure the battery connector and charging port are dry.

Improper wearing may cause the device to malfunction or affect balance, especially during strenuous outdoor activities, increasing the risk of falls and injury.

Before use, make sure the battery clicks firmly into place. An incorrectly installed battery may loosen or fall out during operation,creating a risk of falls or injury.

Before inserting the charging cable, ensure that the battery connector and charging port are dry.

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada' s licence-exempt RSS(s). Operation is subject to the following two conditions:
1.This device may not cause interference.
2.This device must accept any interference, including interference that can cause undesired operation of the device.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:
1. This device may not cause harmful interference.
2. This device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

⚠WARNING!

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:
-- Reorient or relocate the receiving antenna.
-- Increase the separation between the equipment and receiver.
-- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
-- Consult the dealer or an experienced radio/TV technician for help.

FIRE & ELECTRIC SHOCK HAZARDS

1. Only use the battery pack and charger supplied with this product. Use of any other battery or charger may cause fire, explosion, leakage, or overheating.
2. Do not charge or use the product if the power cord, plug, or battery is damaged, frayed, wet, swollen, leaking, or emits an unusual odor or heat.
3. Do not use the product in the presence of flammable or explosive gases, vapors, or liquids.
4. This product contains a rechargeable lithium-ion battery. Do not dispose of it in fire or with household waste. Recycle according to local regulations.

PERSONAL INJURY HAZARDS

5. It is not a medical device and does not prevent falls or replace medical treatment.
6. This product is not intended for use by children. Keep the product, battery, and charger out of reach of children and pets.
7. Do not exceed the maximum user weight or usage limits specified in the technical parameters.
8. Always ensure all straps are securely fastened and properly adjusted before operation. A loose fit may cause loss of balance or equipment damage.

⚠WARNING!

9. Do not submerge the product in water or use it in heavy rain beyond its stated IP rating. Do not wear it while bathing or swimming.
10. Do not modify, disassemble, or repair the product yourself. Repairs must be performed by authorized service personnel.
11. Inspect the product before each use. Do not use if any part is cracked, deformed, or malfunctioning.
12. Use only manufacturer-approved accessories and replacement parts.

Additional Safety Warnings

Read these warnings and all other safety instructions before use.

1. Users who must not use this product

Do not use this product if you have difficulty walking or standing, or cannot maintain your balance independently.
Do not operate or use it after drinking alcohol or taking medication that impairs judgment or reaction time.

2. Prohibited environments and activities

Do not use the product near strong electrostatic discharges or strong magnetic fields, which may cause electronic system malfunctions.
Do not use it for high-risk sports or on complex or hazardous terrain where secure footing or balance cannot be maintained, including unstable slopes, slippery surfaces or unprotected edges.
Outdoors, assess the terrain, weather and your physical condition, stay within your abilities, and take responsibility for your safety. This product does not prevent falls.

3. Unexpected shutdown or abnormal operation

During normal use, stop if the device stops unexpectedly, produces unusual sounds, abnormal high-pitched noise or vibration, or displays signals inconsistent with the normal indicators described in this manual. Stop moving as soon as it is safe, establish a stable position, switch off the device and remove the battery if safe. Do not restart or continue use; contact authorized technical support and have the device checked by authorized service personnel before reuse.

⚠WARNING!

Additional Safety Warnings (continued)

4. Risks of incorrect wearing

Follow this manual: align the motors with your hip joints, securely fasten the waist buckle and leg straps over clothing, and maintain a comfortable fit without excessive tightness. Activate assistance only after correct fitting and successful wear detection.

Misalignment, loose or unsecured straps or excessive tightness may cause skin irritation or pinching, impaired circulation, restricted movement, excessive forces on the hips or legs, loss of balance, falls and injury. If you experience pain, numbness, slipping or restricted movement, stop use and switch off before adjusting or removing the product.

5. Risks of incorrect battery installation or charging

Switch off before installing or removing the battery. Align its guide slot with the rail at the left front mounting point; slide it in until the latch clicks and check that it is locked. Do not force it. Incorrect installation may cause detachment, power interruption, unexpected loss of assistance, falls or injury, and damaged electrical contacts.

Remove the battery for charging. Use only the supplied battery, charging dock, charger and cable. Keep connectors and the charging port clean and dry. Charge at 3°C to 35°C in a dry, ventilated area away from flammable materials. Let a hot battery cool first. Unplug the charger and cable when charging is complete.

Do not charge a battery that is damaged, wet, swollen, leaking, unusually hot or emitting an unusual odor, or use charging equipment in any such condition. Incorrect installation or charging, incompatible equipment, moisture or short circuits may cause overheating, leakage, fire, explosion, electric shock, burns, injury or property damage.

Users are responsible for following the wearing, battery installation and charging instructions and for safe operation and use. These warnings do not exclude or limit rights or liabilities that cannot be excluded or limited under applicable law.

Legal Liability and Disclaimer

The terms and conditions described in this chapter apply to Reudyn's wearable exoskeleton products. Models or versions not explicitly listed herein are not covered by this statement; applicable provisions may be found in other documents.

Reudyn is intended solely for outdoor recreational use (specific usage restrictions are detailed in the Safety Guidelines). It is not designed for medical, rehabilitation, or diagnostic purposes, nor can it be used to diagnose, alleviate, treat, or prevent any disease or medical condition. Reudyn shall not be liable for any injuries, damages, or other consequences resulting from use beyond the stated purposes or improper operation.

By using this product, you acknowledge that you have read, understood, and accepted all terms and conditions contained herein, and agree to strictly comply with all product usage specifications. Failure to adhere to these provisions—including those outlined in this statement, manual, safety guidelines, or any related documentation—will result in consequences for which the user shall be solely responsible.

For personal injuries or property losses directly or indirectly caused by improper use, unauthorized modification, failure to heed warnings, or similar reasons, the user shall bear full responsibility.

To the greatest extent permitted by law, Reudyn's liability shall be limited solely to repairing or replacing products with confirmed material, workmanship, or design defects in accordance with Reudyn's valid warranty policy. Reudyn shall not be liable for any indirect, incidental, consequential, or punitive losses—including but not limited to loss of profits, data loss, business interruption, damage to goodwill, or opportunity losses—regardless of their cause (including product defects) or whether Reudyn was notified of such losses. Reudyn shall not be liable for any consequences arising from the use of unauthorized or unapproved third-party components or unauthorized modifications.

Notwithstanding the aforementioned provisions, nothing in this clause is intended to exclude or limit any liabilities that Reudyn is legally obligated to assume under applicable law and which cannot be restricted by agreement.

Trademarks and Copyrights

Reudyn and related marks are trademarks of Reudyn Innovation (Shenzhen) Technology Co., Ltd. (hereinafter "Reudyn") and its affiliated companies. All other trademarks, service marks, and trade names mentioned in this manual belong to their respective rights holders.
Copyright © 2026 Reudyn Innovation (Shenzhen) Technology Co., Ltd.
All rights reserved. Without prior written permission from Reudyn, no part or all of this manual, its accompanying documents, or software may be reproduced, stored in retrieval systems, translated, or distributed in any form or by any means (including but not limited to electronic, mechanical, photocopying, audio recording, or other methods).
Information in this manual is subject to change without prior notice. Reudyn reserves the right to revise or modify this manual at any time.

Reading Guidelines

This manual is intended for users and personnel responsible for operating and maintaining the product. Read it carefully before use and keep it for future reference.

⚠ CAUTION

Indicates a potential hazard that, if not avoided, may cause minor personal injury or equipment damage.

ⓘ Notice

Highlights important information or practical tips that do not involve a risk of personal injury.

Download the Reudyn App


You can download the latest version of the Reudyn app using one of the following two methods.

- Method 1: Click the link below to download the app.


 https://reudyn.com/pages/download-center

- Method 2: Search for "Reudyn" in the App Store or App Store China to download the app.


WHAT'S IN THE BOX

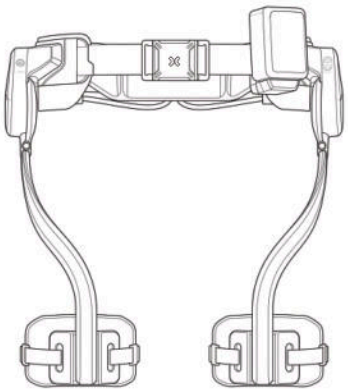Main Unit × 1

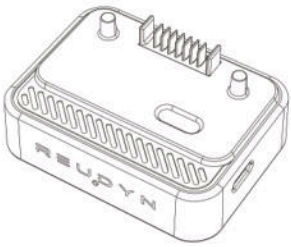
Charging Dock × 1

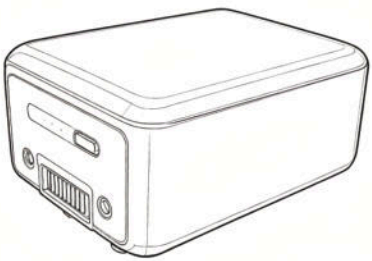Battery ×1


Manual ×1

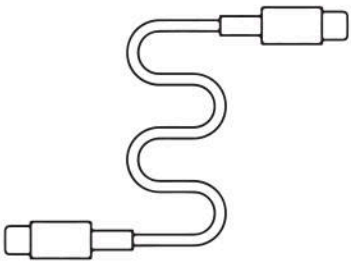Charging cable ×1

PRODUCT OVERVIEW

Diagram

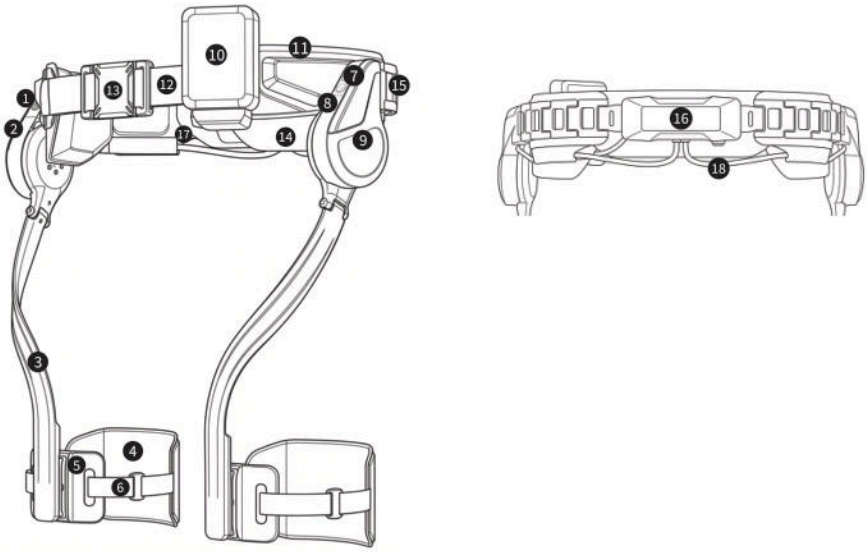1	Power Button	7	Mode Button	13	Waist Buckle
2	Status/Battery Indicator	8	Gear Indicator Lights	14	Back Strap
3	Leg Bar	9	Charging Port	15	Waist Strap
4	Rear Leg Strap	10	Battery Compartment	16	Main Module
5	Leg Band	11	Waist Pad	17	Battery Release Button
6	Leg Buckle	12	Waist strap	18	Belt Adjustment Button

The image is for reference only. Actual product design, features, and specifications may change without notice.

Specifications

Power Form: Electric Energy Storage
Folded Dimensions: 390 × 295 × 170 mm
Weight: About 2.6kg
Battery Weight: 400 ± 30 g
Rated Voltage: 21.6 V⎓
Rated Capacity: 2.8 Ah
Nominal Energy: 60.48 Wh
Peak Output Power: 1000 W
Peak Torque: 2 × 21 N·m
Range: 20 km
Recommended User Height: 165–185 cm
Waist Size Compatibility:
Size M: 80–100 cm
Size L: 95–120 cm
Recommended Leg Circumference: 40-62 cm (below the femur and above the knee joint)
Water Resistance: IP54
Operating Temperature: −10°C to 40°C
Materials: Carbon Fiber, Magnesium Alloy, Rubber
Designed For: Hiking, Climbing, Cycling, Photographers, Training and more
Operating Frequency: 2402 - 2480 MHz
Maximum transmission power: 2.67 dBm
Channel Spacing: Low Energy mode: 1MHz, 2MHz

WEARING AND REMOVAL


Wearable Device

Proper fitting of the exoskeleton—from the waist down to the legs—ensures optimal performance without restricting mobility.


First, rotate both leg bars outward to unfold the device.

Attaching the Waist Belt

1.Unlock the left and right waist straps. Pull the belt adjustment button to one side; both straps will extend simultaneously. Release the button to lock them automatically.

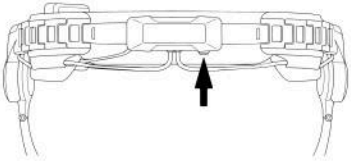
2.Position the belt so the main module rests at the lumbosacral junction (level with the top of the hip bones/iliac crest). When viewed from the front, the waist buckle should align approximately with the navel. The hip joint motor should align with your hip joint.

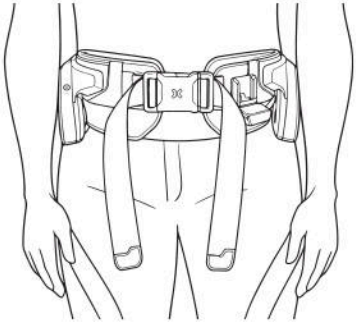3.Firmly fasten the buckle until you hear a click. Cross the straps behind you: use your right hand to pull the left strap backward and your left hand to pull the right strap backward, securing the belt snugly so it does not shift during movement.

4.When properly worn, the motor should be aligned with the hip joint and maintain a gap of 0.5–1 cm.

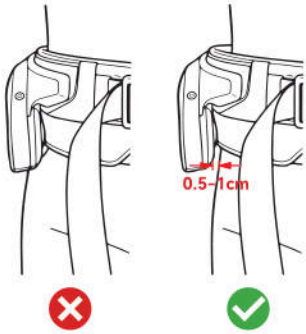
⚠CAUTION

• Ensure both upper and lower snap tabs are fully inserted and securely fastened.


• When fastening the waist buckle, ensure it does not pinch clothing.

ⓘNOTICE

• If the motor sits too close or too far from the hip bone, adjust the lumbar bar width until the motor aligns with the hip joint and maintains a 0.5–1 cm gap.

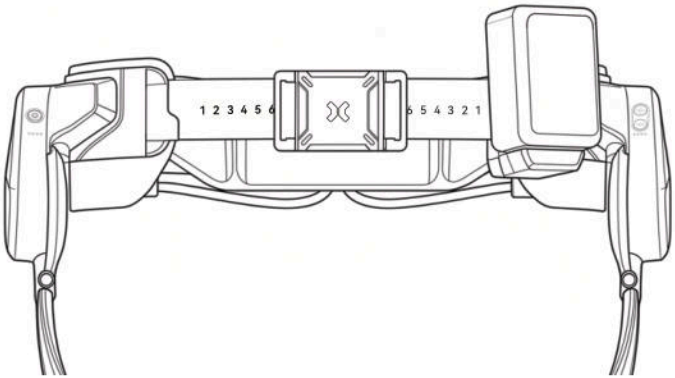
• If the strap is difficult to pull, use your other hand to hold the belt in place before tightening


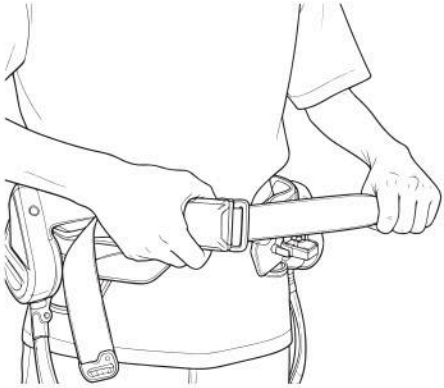
Attaching Leg Straps

1.Ensure the device is centered on your body's central axis for balanced performance.

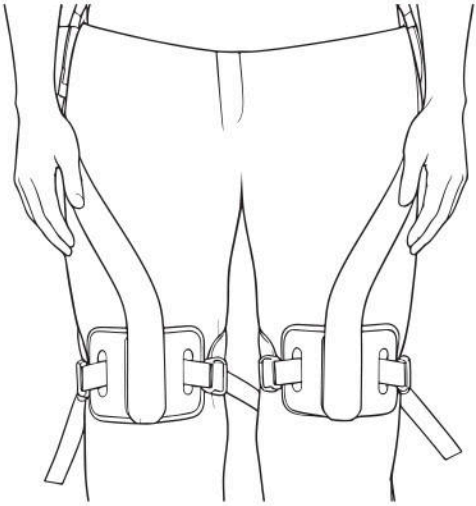
2.Position the rear leg strap behind your thigh, just above the knee. Ensure the rear leg strap does not rest in the popliteal fossa (back of the knee) or interfere with knee flexion and extension

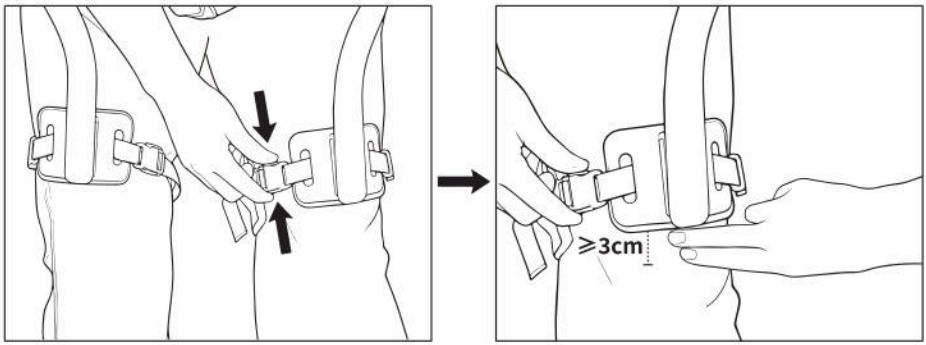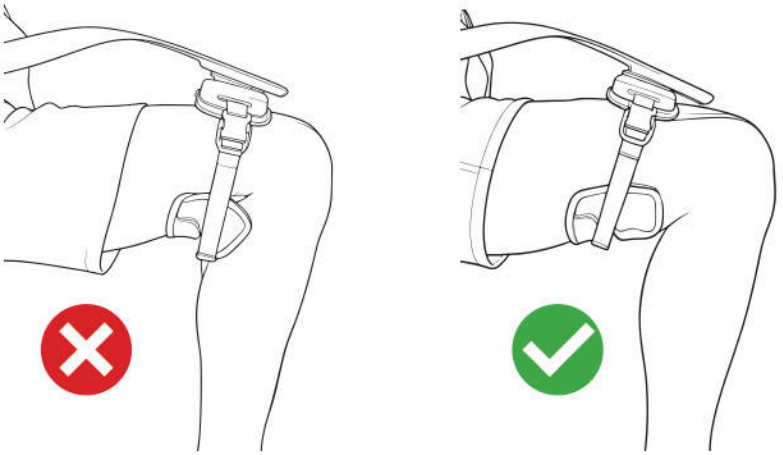3. The leg band should sit horizontally centered on the back of the thigh. Adjust the leg straps for a comfortable fit—both left and right segments on each leg are independently adjustable.

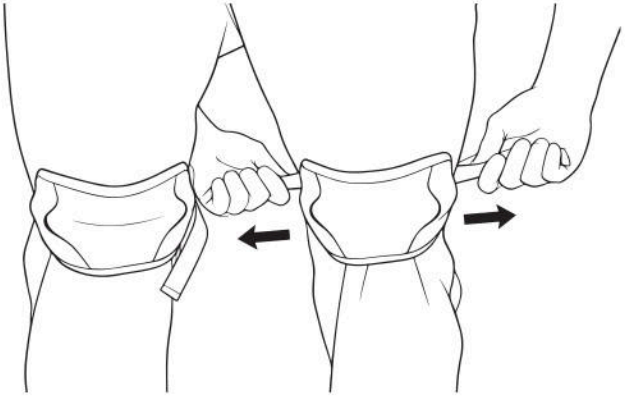4. Ensure all straps are correctly positioned and fully tightened.

⚠CAUTION

• Tighten the strap securely while maintaining comfort; excessive tightness may restrict blood circulation. Ensure the leg band fits snugly against the thigh without compressing tissue or causing discomfort.
• During activity, ensure the waist belt and leg straps are worn over clothing to prevent injury. Avoid direct contact with exposed skin, as this may cause irritation or allergic reactions.
• Avoid wearing skirts or excessively loose clothing to ensure proper fit and function.

ⓘNotice

• When wearing long pants, it is recommended to pull pant legs slightly upward to prevent the fabric from catching.
• During walking, the waist belt and leg straps may be worn moderately snug; during high-intensity activities such as running or cycling, all straps must be fully tightened for stability and optimal support.

Removing the Device


⚠CAUTION

To avoid safety hazards, please power off the device before removal.


ⓘNotice

To prevent the exoskeleton from falling accidentally when removed, it is recommended to hold the back frame with one hand while unfastening the waist buckle with the other.


1. Press and hold the control button to shut down the device.

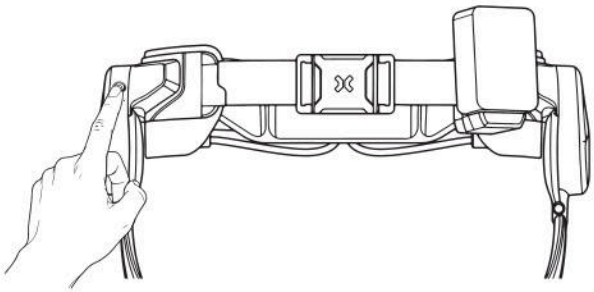2. Press and hold the buttons on both sides of the leg buckles to remove the leg straps.

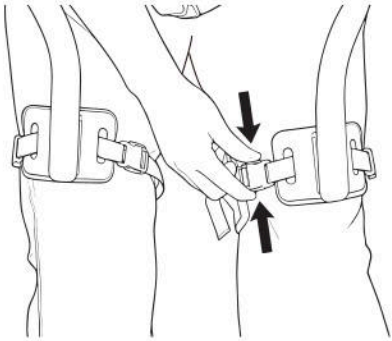
3.Hold the waist button, press it toward the center, and then pull the top spring tab upward.

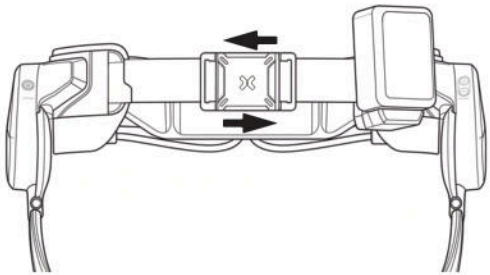
Stowing the Device

1.Press the waist adjustment buttons on both left and right lumbar bars. Push both lumbar bars inward to their minimum width setting.

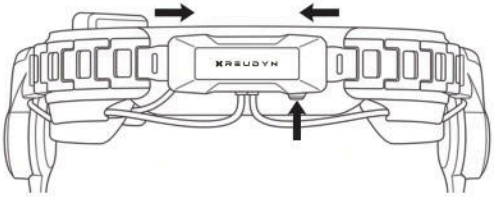2.First, lift the left leg support and fold it inward. Then, lift the right leg support and fold it inward.

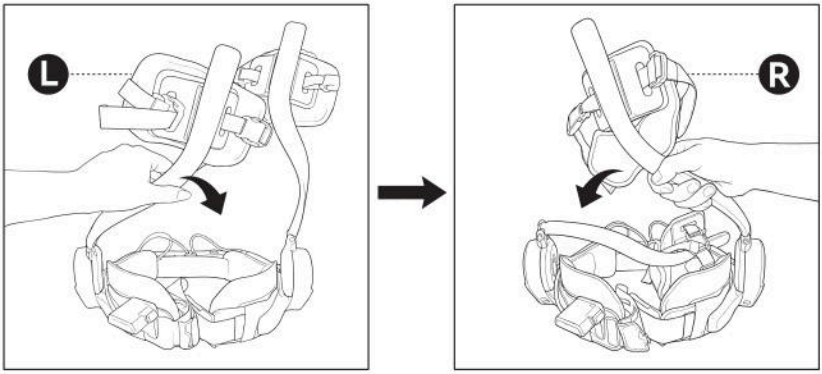3. Grasp the intersection point of both leg bars, rotate outward, and place the device into the carrying case for storage.


FIRST USE


First-Time Charging

For safety, the battery is shipped in sleep mode. Before first use, it is recommended to fully charge the battery for optimal performance.

1. Remove the battery from the box; it is packaged separately.

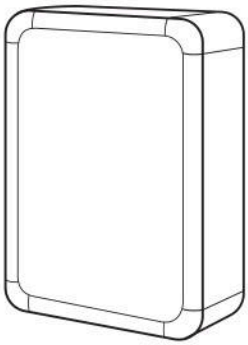2.Connect the battery to the charging port of the included charger and connect the adapter to the power source.

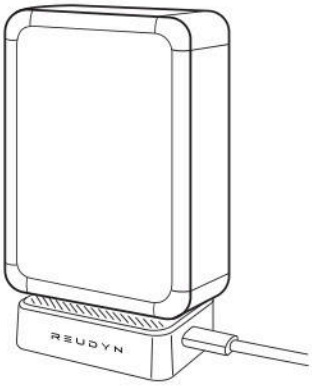⚠CAUTION

• For use only with C1387B6 battery (battery pack, etc.)

Power On/Off

Power on

Long-press the right control button once; you will hear a beep and both status lights will flash green. The device is now powered on.

Power off

Press and hold the right control button once to power off the device.

ⓘNotice

Before powering on the device, ensure you are wearing it correctly.

Connecting the Device

1. Open the Reudyn app and log in using your phone number and verification code. If you do not have an account, the verification code will log you in and create a new account automatically.

2.Follow the on-screen instructions to enter your body data. The app will calculate recommended assistance settings based on your height and weight.

3.Long-press the right control button to power on the device. Place your phone near the device and tap Search Device to begin connecting.

4.After the device connects successfully to your phone, automatic wear detection will run, and results will appear in a pop-up window. If detection fails, follow the on-screen instructions to adjust the fit, then tap Redetect to try again. You may only proceed after successful detection.

ⓘ Notice

• The device runs automatic wear detection each time it powers on.
• If you are seated, detection may fail. Stand up and tap to retest.
• Once connected, this exoskeleton will appear in your device list. View it under My Devices (Settings icon in the top-right corner).
• Each device can only connect to one phone at a time. To use a different phone, first disconnect the device from the original phone in Settings > My Devices. Only then can another phone connect.

ASSISTANCE MODES

Smart Mode

AI adaptive mode: The device automatically identifies terrain, activity type, and posture, then adjusts assistance strategy autonomously; users cannot manually adjust settings.

Standard mode

In this mode, users can adjust the exoskeleton's assistance level according to their needs; no automatic scene detection is used.

Custom Mode

Users can use preset assistance profiles in the app or customize parameters to their individual preferences.

OPERATING THE DEVICE

Control Buttons

Power On/Off

Long press for 2 seconds
Long-press the right control button once to toggle the device on or off.

Switching Assistance Modes

Double-press
Press the right control button twice quickly to cycle through modes. Factory default is Smart Mode; cycle order is: Smart Mode → Standard Mode → Custom Mode.
Single-press
Press the right control button once to activate or pause power assistance in the current mode.

Adjusting Assistance Level

Single-press
Press the left adjustment button (+/-) to increase or decrease assistance intensity by 25%.

Rerunning Wear Detection

Single-press
If all status indicators flash orange, check that the device is properly worn, then press the control button once to rerun automatic wear detection.

Status Indicators

The status indicators are located on the left and right motors, arranged from left to right. Adjacent to the control buttons are LED 1, followed by LED 2, LED 3, and LED 4 from left to right.


During Operation

Button	Function	Indicator Behavior
Right Button (Mode)	Mode Switch	Double-click the right control button to cycle through assistance modes in the following order: Smart Mode → Standard Mode → Custom Mode.
		Double-click to cycle modes: • Smart Mode: AI fully adaptive, gear locked at level 4, green light.
		• Standard Mode: Constant power assist, level-adjustable, blue light.
		• Custom Mode: Uses app preset configuration, orange light. If unconfigured at first use, defaults to Fitness Mode.
Left Button (Gear)	Power Adjustment	• Standard/Custom Mode: Press + to upshift, press – to downshift. 4 gear levels available (25% to 100% assistance, in 25% increments). Indicator lights update accordingly.
		• Smart Mode: Pressing +/- does not change gears; all four gear indicators remain solid.

Battery


Removing the Battery

Press and hold the battery release button, then carefully slide the battery out of the waist compartment.

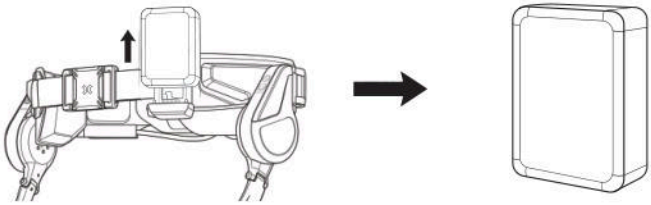
Installing the Battery

Align the guide slot on the battery with the rail on the left front mounting point. Push the battery into the latch until you hear a click, indicating it is securely installed.

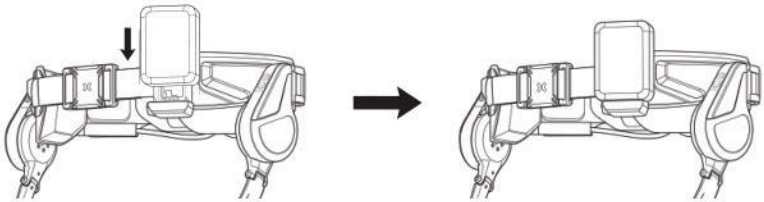Charging the Battery

1. The battery must be removed from the device.
2. Connect the charger to a power source via the charging port. A flashing red indicator light means charging is in progress; a green light indicates a full charge.

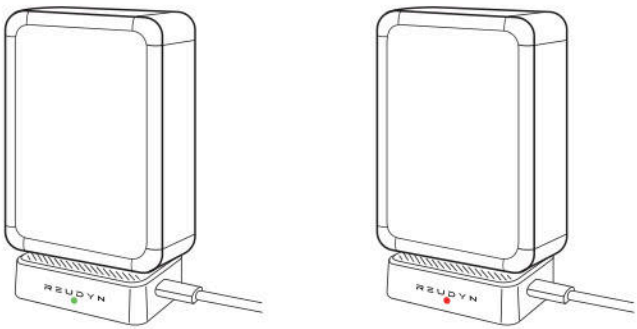3. After the battery is fully charged, unplug the charger and charging cable.

Low-Battery Safety System

The exoskeleton features a built-in low-power safety system designed to ensure your safety.
When battery level drops to a critical threshold, the device will emit an audible alarm to remind you to charge or replace the battery. As battery level continues to drop, the device will gradually reduce assistance. In low-power mode, acceleration mode is unavailable; when battery level drops further, assistance will cease entirely.
As battery level continues to drop, the device will gradually reduce assistance (to no more than 10%) until it runs out of power and shuts down automatically. In low-power mode, acceleration mode is unavailable; when battery level drops to 5%, extreme mode is unavailable.

Reudyn APP

Device status

No.	Description	No.	Description
1	Currently connected device	5	Pause/Resume Assistance
2	Status bar: connection status, battery percentage, distance	6	Current mode & assistance level
3	Settings: voice prompts, device settings	7	Activity type selection
4	Notification Center: system messages	8	Terrain type selection

Firmware Update

Update the firmware to access the latest features.

1. Ensure your mobile device is connected to the internet and make sure the exoskeleton is connected to your phone.
2. On the profile page, go to My Devices > Settings > Check for Updates. If new firmware is available, tap Upgrade to begin the update.
3. If the update succeeds, restart the device. If the update fails, follow the on-screen instructions to retry or contact customer support.

ⓘ Notice

• During the upgrade, ensure your mobile device is connected to the internet stably. Stay in the Reudyn app and do not exit or switch to another application.

• During the upgrade, ensure the battery level is above 20%. It is recommended to perform both the upgrade and charging simultaneously.

• During the upgrade, the device may restart multiple times. Keep the device facing upward; the leg rod may move—take safety precautions.

MAINTENANCE AND STORAGE

Cleaning the Main Unit

Cleaner Unit Structure

• Gently wipe the device with a soft, damp cloth.

• Ensure the exoskeleton is completely dry before storage or use

⚠ CAUTION

• Do not use corrosive or volatile chemical solvents (such as alcohol, gasoline, diesel, acetone) or abrasive materials, as these may damage the surface and internal components.

• Do not use a high-pressure washer or high-pressure water gun on the main unit; avoid exposing the device to pressurized or high-velocity water streams.

• Do not submerge the device in water or allow electronic components to get wet. When cleaning, prevent water from entering the charging port or battery interface.

ⓘ Notice

The lumbar support and adjustment buckle feature a tightly fitted design. Over time, their surfaces may develop light marks from friction. Wipe these marks with a soft cloth dipped in 75% alcohol and allow them to dry naturally.

Cleaning Soft Fabrics

• Wipe soft fabric surfaces with a soft damp cloth and mild detergent.

• Ensure soft fabrics are completely dry before storage or use

⚠ CAUTION

• Do not use detergents containing bleach or strong alkaline substances, as these may damage soft fabrics.

• Do not machine wash, as vigorous agitation may deform or damage the soft fabric.
• Do not tumble dry, and avoid exposing soft fabrics to direct sunlight or high temperatures, as this may cause fading or deformation.

ⓘ Notice

Soft fabrics are consumable parts and will naturally wear during normal use. If wear, tearing, damage, or deformation occurs, replace them promptly.

Battery Maintenance

• Before each use, inspect the battery for damage, bulging, or leakage.
• Charge the battery when its level drops to 15%; avoid overcharging.
• Recommended charging temperature range: 3°C to 35°C. Charging will automatically pause outside this range.
• Do not charge immediately after use; allow the battery to cool down first.
• If the battery will not be used for more than 10 days, recharge periodically to maintain battery level between 40% and 60%.

Storage and Transportation

Storage

• Store the device in a clean, dry, well-ventilated indoor environment; recommended relative humidity is 45–75%.
• Do not store the device outdoors in direct sunlight, rain, dusty, or humid environments.
• Keep away from extreme temperatures, open flames, and heat sources to prevent material degradation.
• Do not place heavy objects on the device to prevent deformation or mechanical damage.
• For long-term storage, cover the device with a protective cover to prevent dust accumulation or impact

Transportation

• Before transport, ensure the device is properly folded and securely placed in the carrying case.
• Battery transport is subject to applicable aviation regulations and airline policies. Check with your airline in advance for specific requirements.

• During transport, protect the device and accessories from impacts, vibration, or excessive compression.
• During transport, do not expose the device or battery to extreme temperatures.
• When transporting over long distances, ensure the device is completely powered off.
• If packaging shows signs of damage, inspect the device and battery before use to ensure safety.

This manual is subject to change without prior notice. If you have any questions about this manual, please contact us by email:
support@reudyn.com

Copyright © Reudyn. All rights reserved.
Follow us:
Warranty Policy
https://reudyn.com/pages/download-center

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Copyright © Reudyn. Tous droits réservés.